Otovent® GO

EN

Otovent® GO is a medical device that helps you equalise ear pressure easily and conveniently to relieve symptoms of ear pressure, e.g. after a cold, when flying or when diving. Otovent® GO is easy to use and bring, and provides relief whenever symptoms arise.

Contents


A. Handle (polypropylene)

B. Nosepiece, 3 different resistance levels;

01. LIGHT

02. MEDIUM

03. STRONG

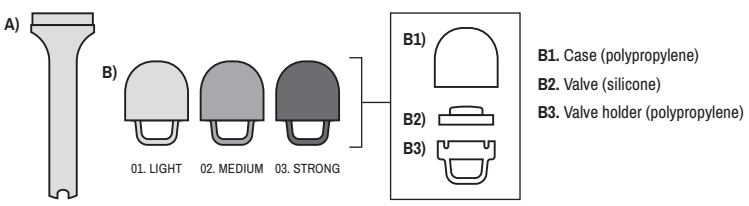B1. Case (polypropylene)

B2. Valve (silicone)

B3. Valve holder (polypropylene)

Otovent® GO is supplied pre-assembled and ready for use.

Instructions for Use

Read the instructions before use.

Indications and clinical benefits

Otovent® GO helps to equalise negative pressure in the middle ear.

Common symptoms of negative pressure in the middle ear:

- Ears feel blocked
- Impaired hearing
- Ear discomfort or pain
- Sensation of heaviness or light dizziness

How does Otovent® GO work?

Otovent® GO helps to equalise ear pressure in the middle ear when air is exhaled through the nose against resistance. This increases pressure in the nasal cavity, which allows the Eustachian tube to open.

Before using Otovent® GO

Remove the Otovent® GO from its packaging by pulling the round end of the nosepiece upwards.

Otovent® GO is supplied pre-assembled and ready for use.

How to use Otovent® GO

General use of Otovent® GO

Otovent® GO should be used three times a day (morning, afternoon, and evening) or more often if needed, with one exhalation per nostril each time. If you do not experience any relief from your symptoms within one week, you should contact your doctor.

For symptom relief during flights, diving, or treatment in a hyperbaric chamber

When flying, use Otovent® GO 3–4 times as soon as the aircraft begins its descent or if you experience discomfort during the flight. If your ears continue to feel blocked, you can continue to use Otovent® GO after landing. Otovent® GO can be used multiple times to equalise ear pressure before/after diving or during/after treatment in a hyperbaric chamber.

User Instructions

Otovent® GO is supplied pre-assembled and ready for use, equipped with the **“01. LIGHT”** nosepiece (low resistance).

1. Position the round part of the device with the small round opening firmly against one nostril.
2. Use your index finger to close your other nostril.
3. Inhale deeply through your mouth, then close your mouth and exhale through the Otovent® GO until resistance reduces and air passes through the device.
4. Repeat the procedure through your other nostril.

Successful use is often characterised by a distinct “pop” or the feeling that ears have been unblocked.

If no effect is achieved with resistance **“01. LIGHT”**, pull apart nosepiece and handle and switch to the nosepiece marked **“02. MEDIUM”** and repeat the procedure. If this still does not provide the desired result, switch to the nosepiece marked **“03. STRONG”** and try again.

It may take several attempts before you master the technique.


Cleaning

1. Remove the handle by pulling the handle and nosepiece apart.
2. Separate the valve holder from the case (but leave the valve in place).
3. Rinse the case with water. If necessary, rinse the valve holder as well. Ensure that the valve holder is held with the handle facing down while rinsing to prevent the valve from coming loose.

4. Wipe the inside of the case dry.
5. Leave to air dry completely before reassembling.

Disassembly/Assembly

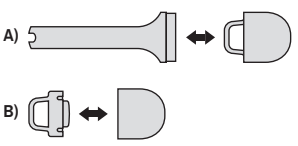Clean only one nosepiece at a time to avoid mixing them up.

If the valve comes loose from the valve holder, reattach it.

Contraindications

When should Otovent® GO not be used?

Otovent® GO should not be used in cases of acute ear infection (acute otitis media), upper respiratory tract infections, if you have a history of eardrum problems, or if you have recently undergone sinus surgery.

Precautions

If you suffer from frequent minor nosebleeds, Otovent® GO should be used with caution since any manipulation of the nose may cause a nosebleed.

Consult your doctor before using Otovent® GO if you:

Have had a perforated eardrum, suffered facial trauma/fracture, or undergone surgery involving the ears, facial bones or sinuses.

Suffer from poorly controlled, severe, or unstable glaucoma.

Have a serious cardiovascular condition such as heart failure or arrhythmia, due to the risk of stroke or cardiovascular instability.

Stop using this product and contact a doctor if you experience:

Severe headache, severe dizziness, nausea, sensivity to light or severe discomfort from the ear.

Side effects

Mild and temporary side effects may occur, such as temporary discomfort or pain in the ear, headache, dizziness, or nosebleeds (especially if you are prone to nosebleeds or have recently blown your nose forcefully; this condition is usually harmless).

Additional information

Otovent® GO is suitable for adults and children who are 12 years of age and older.

Children must be treated under adult supervision.

WARNING: Contains small parts - choking hazard. Keep away from children 0-3 years old.

Do not use the product if it is damaged.

If a serious incident occurs during use, contact the manufacturer, info@circiuspharma.com, or the competent authority of your country.

Shelf life, storage, and disposal

Otovent® GO is a reusable product intended for use by a single person only (non-sterile).

Recommended storage: room temperature (18–25 °C) and away from direct sunlight.

Otovent® GO should be used before the expiration date or within 1 year of first use, whichever occurs first.

Dispose of in accordance with local requirements.

Explanation of symbols on the packaging

CE	CE mark
	Manufacturer
	Use-by date (YYYY-MM)
LOT	Batch code
REF	Catalogue number
UDI	Unique device identifier
	Consult instructions for use or consult electronic instructions for use
	Caution
1	Single patient multiple use
0-3	WARNING: Contains small parts - choking hazard. Keep away from children under 3 years old.
MD	Medical Device
	Keep away from sunlight

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Circius Pharma AB